



Wireless Mouse CCTV

User's Guide



Package Contents:


1 each Mattingly Wireless Mouse CCTV

1 each Receiver with Power/video cable

1 each Power adaptor

1 each Carrying pouch

1 each User's Guide

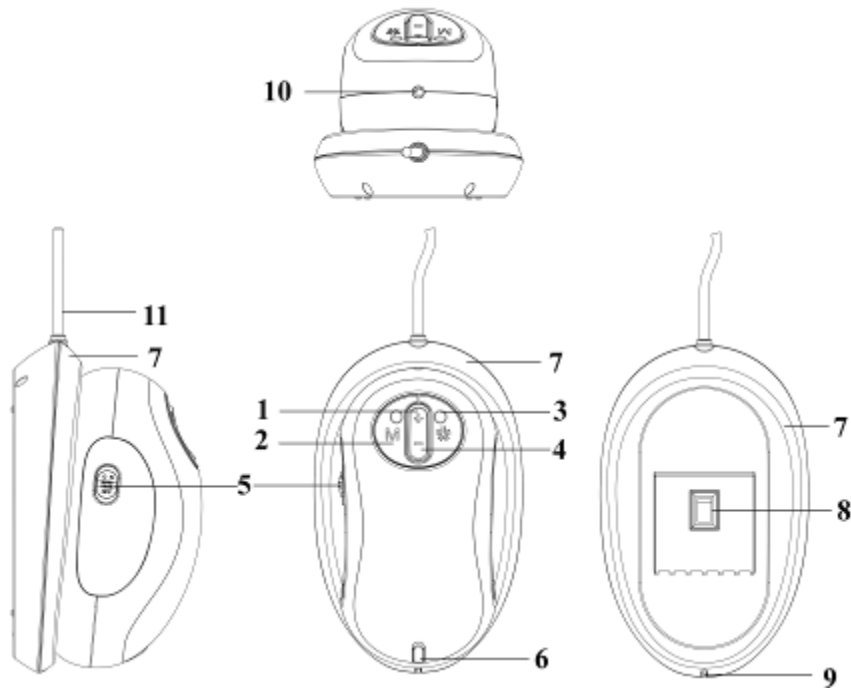

A Quick Look at the Features of the Mattingly Mouse:

1. Increase magnification
2. Mode selection (color, white on black, black on white, or white on black high contrast).
3. Freeze and unfreeze image
4. Lower magnification
5. Mouse power on/off button
6. Mouse charging indicator light (red indicates Mouse is charging; green indicates charging is complete)
7. Receiver (also functions as holder for Mouse when not in use)
8. Receiver power on/off button
9. Receiver power indicator light (red indicates receiver is on)
10. DC-in hole (receives power cable)
11. Y-shaped power/video cable

Charging the Mattingly Wireless Mouse CCTV

The Mouse has a built-in lithium rechargeable battery. Before using the Mouse, charge it until the indicator light turns green (up to 6 hours). If the battery charge is too low, the Mouse will not operate correctly.

1. Turn off the Mouse by sliding its on/off button (Feature#5) to “off”.
2. Plug power adaptor into a wall socket.
3. Plug the end of the power adaptor into the Mouse DC-in hole (Feature#10).
4. The charging light will be red when the Mouse is charging; green indicates the charging is complete (Feature#6).
5. The Mouse will operate for several hours and then will need to be charged again (follow steps#1-4 above). Note: while the Mouse is in use, the charging red/green light will not be on. When the Mouse is ready to be charged again, the picture on the TV will begin to diminish in quality. However, you may recharge the Mouse at any time you wish.
6. Note: the power adaptor can only be used for one purpose at a time: either it is plugged into the Mouse to charge it, or it is plugged in to the Receiver to power it. If you find that having a second power adaptor would be helpful, one can be ordered for you.

How to Use the Mattingly Wireless Mouse CCTV

1. Remove the power adaptor from the Mouse, and plug it into the black end of the Y-shaped power/video cable that is attached to the Receiver (Feature#11). This will provide power to the Receiver. Put the Receiver power button (Feature#8) in the on position; the power indicator on the Receiver will light (Feature#9). Note: the Receiver does not charge the Mouse; it only sends the wireless signal to the Mouse. It can also hold the Mouse when the Mouse is not being used.
2. Insert the yellow end of the Y-shaped power/video cable into one of your television's yellow video-in ports. (If needed, refer to your television owner's manual for more information about locating the video-in port). Note: if your television does not have a video-in port, an alternative RF Modulator can be ordered for you.

3. Turn on your television. If a picture from the Mouse does not appear, switch your television to the appropriate input channel (for example TV “extension” mode, or AV2).
4. Set the Mouse to the magnification level needed (Features#1 and #4), and the Mode needed (Feature#2).
5. Slide the Mouse over the reading material while watching the TV.

Technical Specifications

1. There are eight magnification levels: about 15X, 17X, 19X, 22X, 25X, 30X, 37X, and 50X on a 20 inch TV. The magnification levels are proportional to the size of your television screen. For example, if your TV is twice as large (40 inch) as a 20 inch TV, magnification levels are doubled.
2. There are four modes: color, black on white, white on black, and white on black high contrast.
3. The freeze and unfreeze feature allows the user to capture as a photo whatever the Mouse is viewing. This photo’s level of magnification and mode can then be changed by using the appropriate features on the Mouse.

Safety

1. The Mouse is not waterproof, so please keep it dry.
2. The white LEDs which illuminate the Mouse target field are very bright. Do not stare directly into these LEDs as they may damage your eyes.
3. Use only the Power adaptor that came with the Mouse. Other adaptors may damage the unit.
4. When battery charging is full, disconnect the adaptor from the Mouse.

Warranty Card

Your **Wireless Mouse CCTV** has a one year guarantee against any manufacturing defect!

If you have a defect in your Mouse, please contact the place of purchase or fill in this card and mail it with the Mouse to the address below.

Date of purchase:		Product No.	
Seller's info	Name:		
	Address:		
	Tel:		
Defect description			
Cause			
User's info	Name:		
	Address:		
	Tel:		

Distributed by:



135 Rio West Dr., Bldg. B

El Paso, TX 79932